




EP-550™ / EP-650™

Inversion Table Assembly Instructions

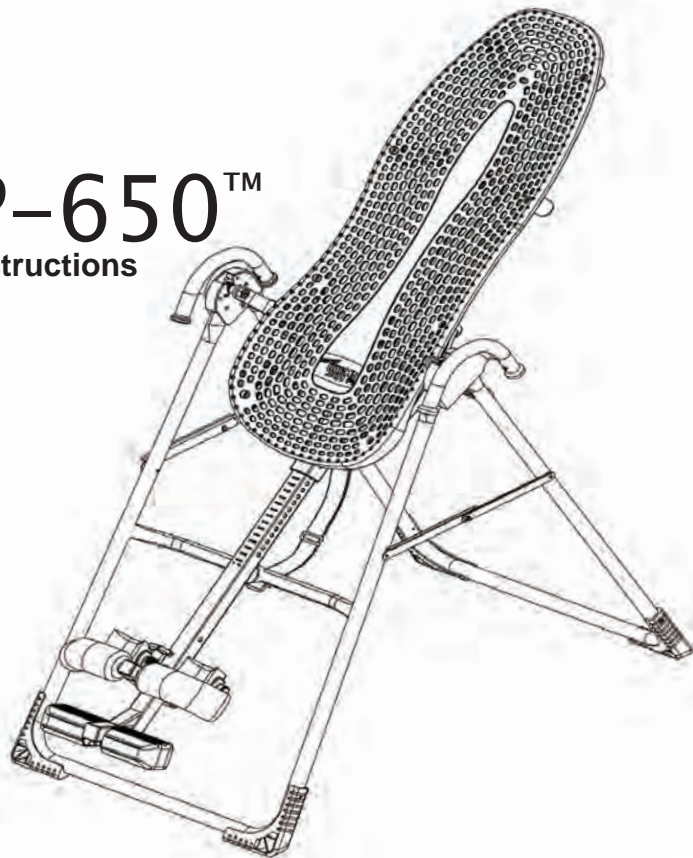

BEFORE YOU BEGIN: Review all steps before beginning assembly and read all precautions before using the inversion table. Carefully adhere to the Assembly Instructions and Owner's Manual to help ensure safety and product integrity.

IMPORTANT SAFETY INSTRUCTIONS
READ ALL INSTRUCTIONS BEFORE USING THE INVERSION TABLE.

 **WARNING**

WARNING - To reduce the risk of injury to persons:

- Read and understand all the instructions, view the instructional DVD, review all other accompanying documents, and inspect the equipment before using the inversion table. It is your responsibility to familiarize yourself with the proper use of this equipment and the inherent risks of inversion, such as falling on your head or neck, pinching, entrapment, or equipment failure. It is the responsibility of the owner to ensure that all users of the product are fully informed about the proper use of the equipment and all safety precautions.
- Close supervision is necessary when the inversion table is used near children, or by or near invalids or disabled persons.
- Use the inversion table only for its intended use as described in this manual. **DO NOT** use attachments not recommended by the manufacturer.
- **NEVER** drop or insert any object into any opening.
- **DO NOT** use or store product outdoors.
- **DO NOT** use if you are over 6 ft 6 in (198 cm) or over 300 lbs. (136 kg). Structural failure could occur or head/neck may impact the floor during inversion.
- **DO NOT** allow children to use this machine.
- Keep children, bystanders, and pets away from machine while in use.
- Keep body parts, hair, loose clothing and jewelry clear of all moving parts.
- The inversion table has no user serviceable parts.
- This product is intended for indoor home use only. **DO NOT** use in any commercial, rental or institutional setting.
- **DO NOT** use the equipment without a licensed physician's approval and a review of the medical contraindications, as noted in the Owner's Manual.
- Failure to assemble and/or use the equipment as directed may void the manufacturer's warranty on this product and could result in injury or death.
- Choose a level surface for assembling and operating the table.
- Follow each step in sequence. **DO NOT** skip ahead.
- Make sure that all fasteners are secure.
- **ALWAYS** test and inspect the table. Make sure the table rotates smoothly to inverted position and back.
- **ALWAYS** replace defective components immediately and/or keep the equipment out of use until repair.

FAILURE TO FOLLOW INSTRUCTIONS AND WARNINGS COULD RESULT IN SERIOUS INJURY OR DEATH.

SAVE THESE INSTRUCTIONS


Items for Assembly

Carefully remove the individual parts from the carton. You should have all of the items listed below. If any items are missing or damaged, contact Customer Service at 1-800-847-0143.


Items for Assembly

Item Numbers


A-frame Base	(F7-1002)
Tether	(F5-1007)
Table Bed Assembly	(EP-1140)
Bolts and Nuts	(EP-1149)
Head Pillow	(EP-1105)

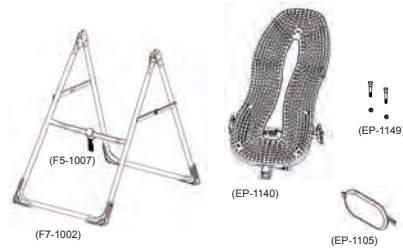

Main Shaft with	
Ankle Lock System	(EP-1120)
Foam Rollers (2)	(F5-1051)
Rubber Plugs (2)	(F5-1056)

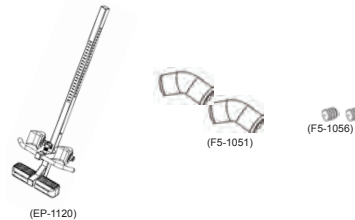


Hand Grips	(F5-1069)
1/2" (13mm)	
hex bolts (6)	(H1-1202)
Roller Hinges (2)	(F5-1064)
Wrenches (2)	(F5-1088)
Allen Wrench	(IA-1149)

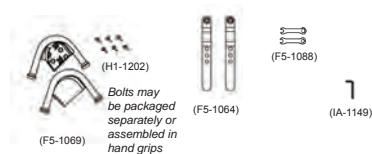

Before Beginning

Before reading further, study the drawing below to familiarize yourself with the important components of your new Teeter Hang Ups® inversion table.

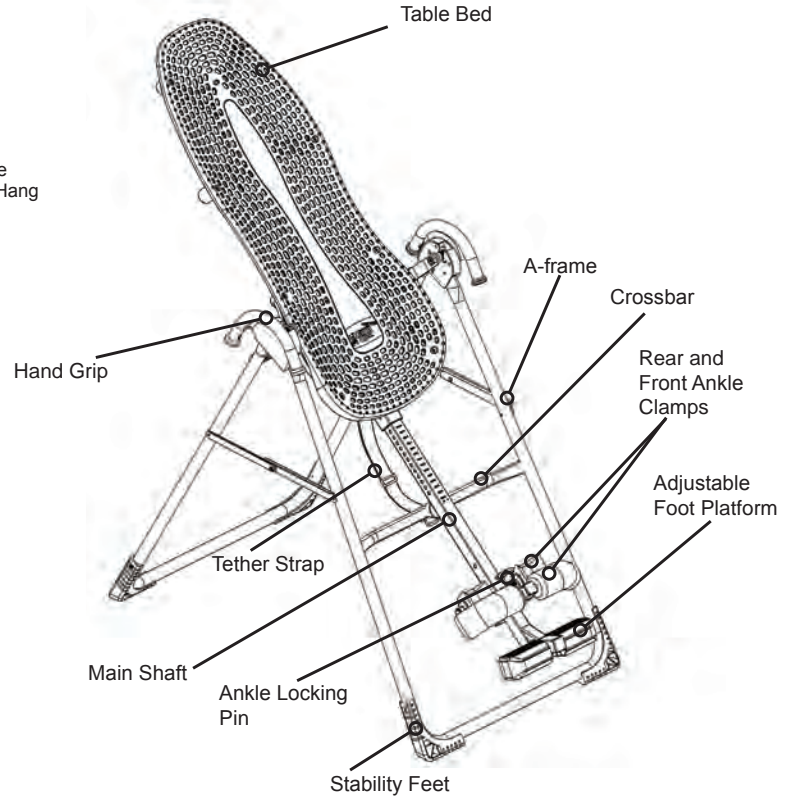

Step 1

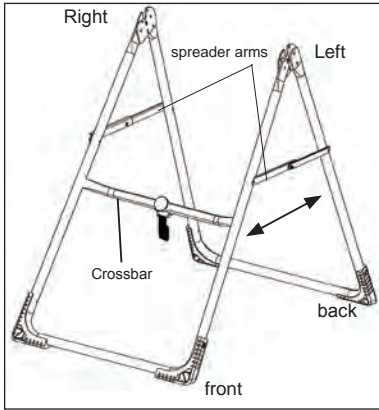


Figure 1

Assemble the A-frame Base (F7-1002) and Install the Handles (F5-1069)

- Open the A-frame and make sure the spreader arms are locked, rest the A-frame on the floor (Figure 1). The Crossbar is at the front.
- Place each Hand Grip (left/right) over the outside edge of the corresponding Hinge Plate (Figure 2).
- Insert three 1/2" (13mm) hex bolts (H1-1202) through the hinge plate into each Hand Grip. Use the wrenches provided to tighten the bolts, being careful not to over-tighten.

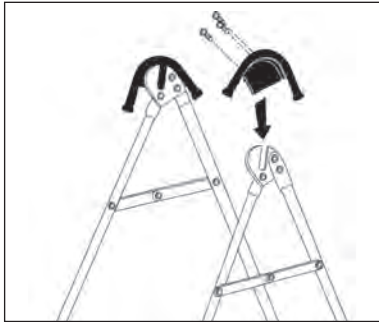


Figure 2

Step 2

Attach the Upper Support Arms to the Table Bed (EP-1140)

- Before following this step, make sure you have read the orange insert that came attached to the Table Bed.
- Place the Table Bed face down on the floor and push down on the Support Beam, so the two holes align evenly with holes at the base of Upper Support Arms (Figure 3). You may have to exert extra pressure to ensure the Support Beam slides over the rubber spacers.
- Insert the two bolts (EP-1149) into the open holes.
- Fasten each bolt with a nut.
- Use the Allen Wrench (IA-1149) to steady the bolts, tighten the nut onto the bolt with the larger wrench (F5-1088) (Figure 4). Repeat with the second bolt.

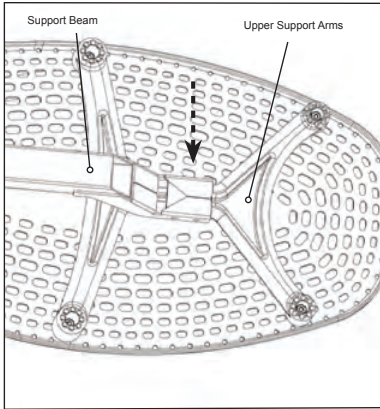


Figure 3

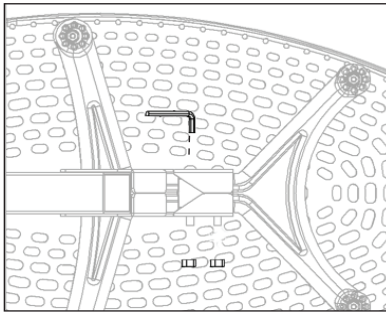


Figure 4

Step 3

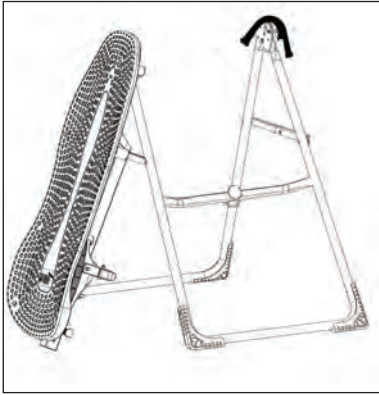


Figure 5

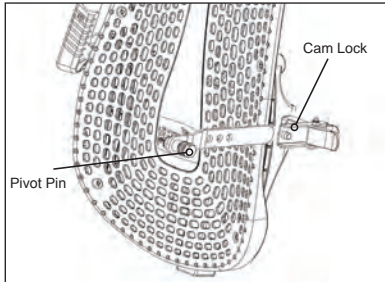


Figure 5A

Install the Roller Hinges (F5-1064) onto the Table Bed (EP-1140)

- For ease-of-assembly, rest the Table Bed against the side of the A-frame (Figure 5).
- Open the Cam Locks on each side of the Table Bed (Figure 5A).
- With the Pivot Pins facing outward, insert the Roller Hinges into the brackets on each side of the Table Bed. The Roller Hinges must slide between the Cam Locks and the brackets.
- Make sure the Roller Hinges are in the same hole setting on both sides. Figure 6 shows the Roller Hinges engaged correctly in setting C.
- Push down on the Cam Locks (Figure 6A) to secure the hinges.

NOTE: Refer to the Owner's Manual for an explanation of the hole settings. If you are unsure, use setting C to start.

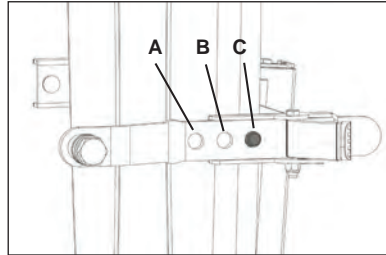


Figure 6

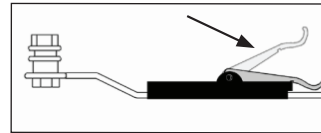


Figure 6A

Step 4

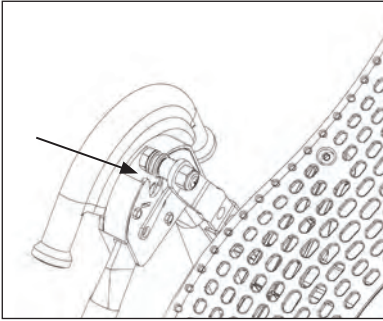


Figure 7

Attach the Table Bed (EP-1140) to the A-frame (F7-1002)

- Holding each side near the Roller Hinges, pick up the Table Bed and stand at the front of the A-frame where the Crossbar is located. Lower each Pivot Pin into the A-frame hinge plates one at a time (Figure 7). You may need to push outward on the A-frame in order to slide the second Pivot Pin into the hinge. Figure 8 shows the correct placement of the Pivot Pin into the hinge plate.
- Make sure that both Pivot Pins are seated at the base of the slot in the hinge plate. Check to make sure that the self-locking hooks have closed over both Pivot Pins, and the table rotates smoothly (Figure 9).

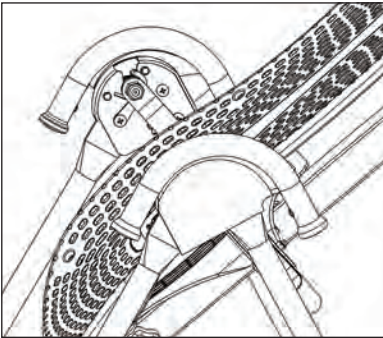


Figure 8

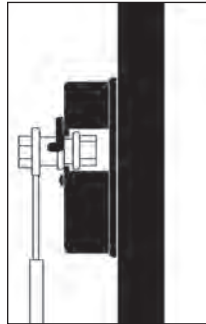


Figure 9

Step 5

Insert the Main Shaft (EP-1120) into the Table Bed (EP-1140)

- Loosen the De-rattler knob on the Main Shaft Housing.
- With the height adjustment settings on the Main Shaft facing up, slide the end of the Main Shaft into the Main Shaft Housing on the Table Bed.
- Pull out the height-selector locking pin to allow the Main Shaft to slide in farther (Figure 10). For the purpose of easy assembly, slide in the Main Shaft and release the pin in the last height setting (Refer to the Owner's Manual for proper height adjustment before use).

NOTE: The De-rattler knob is an optional accessory that helps prevent a slight shift in the table when you invert. The security of the product is not compromised if you choose not to tighten the knob.

- The Main Shaft must rest against the Crossbar of the A-frame (Figure 11). The Crossbar prevents the table from rotating forward when the user steps on the Foot Platform. If the Main Shaft does not rest on the Crossbar as shown in Figure 11 then the Table Bed has been assembled backwards onto the A-frame and this must be corrected before use.

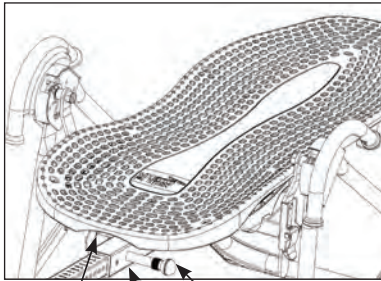


Figure 10

De-rattler Knob

Main Shaft Housing

Height-Selector Locking Pin

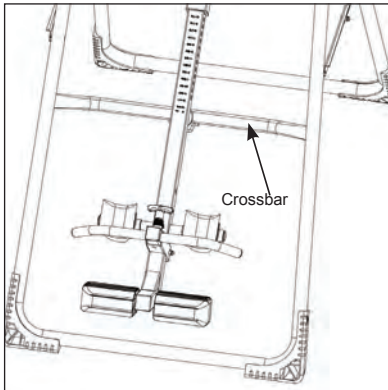


Figure 11

Crossbar

Step 6

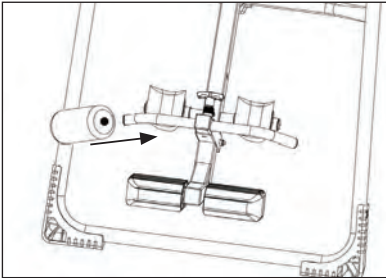


Figure 12

Install the Front Foam Rollers (F5-1051)

- Slide the Foam rollers over each side of the Front Ankle Bar (Figures 12 and 12A).
- Insert the black plugs (F5-1056) into each open end of the Front Ankle Bar if they are not already assembled.

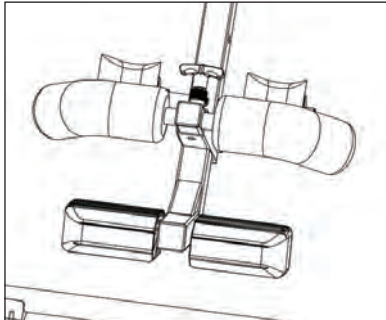


Figure 12A

Step 7

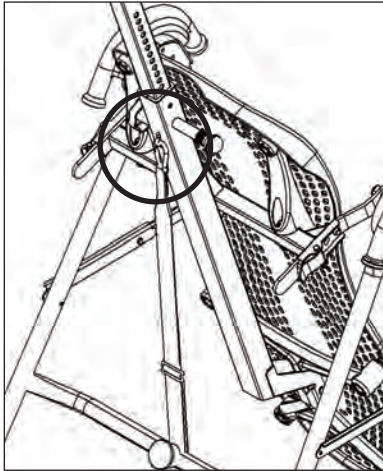


Figure 13

Attach the Tether Strap (F5-1007) to Limit the Degree of Rotation

- Hook the clip at the end of the tether strap to the loop at the base of the table bed (Figure 13) to limit your degree of rotation.
- Tighten or loosen the strap to restrict or extend your angle of rotation.

Attach the Head Pillow

- Attach the Head Pillow by securing the Velcro straps through the holes in the Table Bed – the position of the pillow can be adjusted depending on the user (Figure 14).

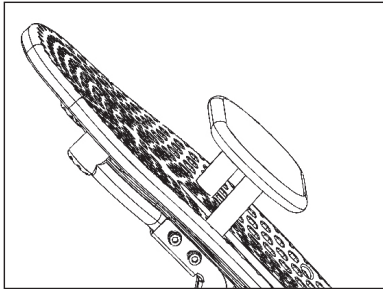


Figure 14

Step 8

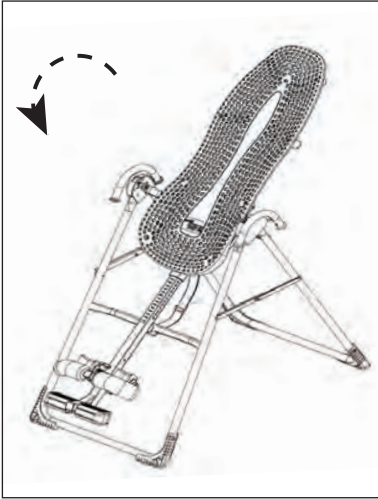


Figure 15

Before Use

- Test the table by hand for smooth and steady rotation (Figure 15).
- Ensure that all fasteners are secure.
- Please complete the warranty registration online at www.teeter-inversion.com.
- For your reference, the serial number can be found on the back of the Table Bed.

WARNING

Read the Owner's Manual thoroughly before using your Teeter Hang Ups Inversion Table.
Improper settings could result in serious injury or death!

Adjustments / Maintenance / Storage



Figure 16

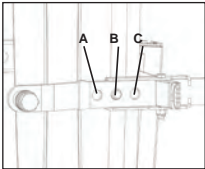


Figure 17

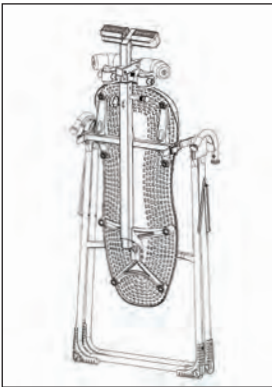


Figure 18

CHANGING THE ROLLER HINGE SETTING

- Stand with your legs on either side of the Main Shaft.
- Reach under each Roller Hinge with your index fingers. Use your thumbs to release the locks over the Roller Hinges (See Figure 16).
- Lift both sides of the table out of the A-frame at the same time. You may rest the table on the Crossbar of the A-frame.
- Unlock the Cam Locks for each Roller Hinge. Change the Roller Hinges to the desired setting (A, B or C) (Figure 17).
- Re-lock the Cam Locks. Replace the Roller Hinges into the hinge plates of the A-frame.

MAINTENANCE

- To clean the Table Bed, wipe down with a damp cloth. Do not use abrasive cleaners.

STORING THE INVERSION TABLE

- Loosen the De-rattler knob.
- Pull the height-selector locking pin and slide the Main Shaft in all the way to the final hole setting. Engage the pin in the storage setting.
- Rotate the table opposite from use until the table has turned 180 degrees and rests against the Crossbar on the A-frame.
- Pull up on the spreader arms to fold the A-frame (Figure 18). **WARNING: This operation may pinch fingers if not done slowly and carefully.**

WARNING

Tipping Hazard: For upright storage, leave A-frame open wide enough to remain stable, or secure to the wall to prevent tipping. In households with small children, the table should be stored flat on the floor, not upright.

CONNECT WITH US!



[facebook.com/teeterhangups](https://www.facebook.com/teeterhangups)



teeter-inversion.com/blog



twitter.com/teeterhangups



[youtube.com/teetertv](https://www.youtube.com/teetertv)

If you have any problems assembling the equipment, or questions about its use, please contact customer service at:

USA: Teeter
Toll Free (Phone) 800-847-0143
(Fax) 800-847-0188
info@teeter-inversion.com
www.teeter-inversion.com

International: Inversion International, Ltd.
(Phone) +1-242-362-1001
(Fax) +1-242-362-1002
info@InversionInternational.com
www.InversionInternational.com



MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany



Any modification to
this device will void
the UL Listing.

U.S. and Foreign Patents Pending. Teeter Hang Ups is a registered trademark of Teeter.
Specifications subject to change without notice.
© COPYRIGHT 2011 Teeter. International Law Prohibits Any Copying

EP-2011_0611-3